240 County Road Ipswich, MA 01938-2723

Tel 978-927-5054 Fax 978-921-1350 www.neb.com info@neb.com

New England Biolabs Product Specification

Product Name: $\alpha 1$ -2 FucosidaseCatalog #:P0724S/LConcentration:20,000 units/ml

Unit Definition: One unit is defined as the amount of enzyme required to cleave > 95% of the α -L-fucose from 1 nmol of Fuc α 1-2Gal β 1-

4Glc-7-amino-4-methyl-coumarin (AMC), in 1 hour at 37°C in a total reaction volume of 10 µl.

Shelf Life: 12 months
Storage Temp: 4°C

Storage Conditions: 50 mM NaCl, 20 mM Tris-HCl, 1 mM EDTA, (pH 7.5 @ 25°C)

Specification Version: PS-P0724S/L v1.0
Effective Date: 11 Mar 2016

Assay Name/Specification (minimum release criteria)

Glycosidase Activity (Endo F1, F2, H) - A 10 μ l reaction in Glyco Buffer 1 containing 1 nM of fluorescently-labeled Endo F1, F2, H substrate (Dansylated invertase high mannose) and 100 units of α 1-2 Fucosidase incubated for 20 hours at 37°C results in no detectable activity as determined by thin layer chromatography.

Glycosidase Activity (Endo F2, F3) - A 10 μ l reaction in Glyco Buffer 1 containing 1 nM of fluorescently-labeled Endo F2, F3 substrate (Dansylated fibrinogen biantennary) and 100 units of α 1-2 Fucosidase incubated for 20 hours at 37°C results in no detectable activity as determined by thin layer chromatography.

Glycosidase Activity (PNGase F) - A 10 μ l reaction in Glyco Buffer 1 containing 1 nM of fluorescently-labeled PNGase F substrate (Fluoresceinated fetuin triantennary) and 100 units of α 1-2 Fucosidase incubated for 20 hours at 37°C results in no detectable activity as determined by thin layer chromatography.

Glycosidase Activity (β -Mannosidase) - A 10 μ l reaction in Glyco Buffer 1 containing 1 nM of fluorescently-labeled β -Mannosidase substrate (Man β 1-4Man β 1-4Man-AMC) and 100 units of α 1-2 Fucosidase incubated for 20 hours at 37°C results in no detectable activity as determined by thin layer chromatography.

Glycosidase Activity (β -N-Acetylgalactosaminidase) - A 10 μ l reaction in Glyco Buffer 1 containing 1 nM of fluorescently-labeled β -N-Acetylgalactosaminidase substrate (GalNAc β 1-4Gal β 1-4Glc-AMC) and 100 units of α 1-2 Fucosidase incubated for 20 hours at 37°C results in no detectable activity as determined by thin layer chromatography.

Glycosidase Activity (β -N-Acetylglucosaminidase) - A 10 μ l reaction in Glyco Buffer 1 containing 1 nM of fluorescently-labeled β -N-Acetylglucosaminidase substrate (GlcNAc β 1-4GlcNAc β 1-4GlcNAc-AMC) and 100 units of α 1-2 Fucosidase incubated for 20 hours at 37°C results in no detectable activity as determined by thin layer chromatography.









New England Biolabs Product Specification

Assay Name/Specification (minimum release criteria)

Glycosidase Activity (β -Xylosidase) - A 10 μ l reaction in Glyco Buffer 1 containing 1 nM of fluorescently-labeled β -Xylosidase substrate (Xyl β 1-4Xyl β 1-4Xyl β 1-4Xyl β 1-4Xyl-AMC) and 100 units of α 1-2 Fucosidase incubated for 20 hours at 37°C results in no detectable activity as determined by thin layer chromatography.

Glycosidase Activity (β 1-3 Galactosidase) - A 10 μ l reaction in Glyco Buffer 1 containing 1 nM of fluorescently-labeled β -Galactosidase substrate (Gal β 1-3GlcNAc β 1-4Gal β 1-4Glc-AMC) and 100 units of α 1-2 Fucosidase incubated for 20 hours at 37°C results in no detectable activity as determined by thin layer chromatography.

Glycosidase Activity (β 1-4 Galactosidase) - A 10 μ l reaction in Glyco Buffer 1 containing 1 nM of fluorescently-labeled β -Galactosidase substrate (Gal β 1-4GlcNAc β 1-3Gal β 1-4Glc -AMC) and 100 units of α 1-2 Fucosidase incubated for 20 hours at 37°C results in no detectable activity as determined by thin layer chromatography.

Glycosidase Activity (α -Glucosidase) - A 10 μ l reaction in Glyco Buffer 1 containing 1 nM of fluorescently-labeled α -Glucosidase substrate (Glc α 1-6Glc α 1-4Glc-AMC) and 100 units of α 1-2 Fucosidase incubated for 20 hours at 37°C results in no detectable activity as determined by thin layer chromatography.

Glycosidase Activity (α -N-Acetylgalactosaminidase) - A 10 μ l reaction in Glyco Buffer 1 containing 1 nM of fluorescently-labeled α -N -Acetylgalactosaminidase substrate (GalNAc α 1-3(Fuc α 1-2)Gal β 1-4Glc-AMC) and 100 units of α 1-2 Fucosidase incubated for 20 hours at 37°C results in no detectable activity as determined by thin layer chromatography.

Glycosidase Activity (α -Neuraminidase) - A 10 μ l reaction in Glyco Buffer 1 containing 1 nM of fluorescently-labeled α -Neuraminidase substrate (Neu5Ac α 2-3Gal β 1-3GlcNAc β 1-3Gal β 1-4Glc-AMC) and 100 units of α 1-2 Fucosidase incubated for 20 hours at 37°C results in no detectable activity as determined by thin layer chromatography.

Glycosidase Activity ($\alpha 1$ -3 Fucosidase) - A 10 μ l reaction in Glyco Buffer 1 containing 1 nM of fluorescently-labeled α -Fucosidase substrate (Fuc $\alpha 1$ -3Gal $\beta 1$ -4GlcNAc $\beta 1$ -3Gal $\beta 1$ -4Glc-AMC) and 100 units of $\alpha 1$ -2 Fucosidase incubated for 20 hours at 37°C results in no detectable activity as determined by thin layer chromatography.

Glycosidase Activity (α 1-3 Galactosidase) - A 10 μ l reaction in Glyco Buffer 1 containing 1 nM of fluorescently-labeled α -Galactosidase substrate (Gal α 1-3Gal β 1-4GlcNAc-AMC) and 100 units of α 1-2 Fucosidase incubated for 20 hours at 37°C results in no detectable activity as determined by thin layer chromatography.


Glycosidase Activity ($\alpha 1$ -3 Mannosidase) - A 10 μ l reaction in Glyco Buffer 1 containing 1 nM of fluorescently-labeled α -Mannosidase substrate (Man $\alpha 1$ -3Man $\beta 1$ -4GlcNAc-AMC) and 100 units of $\alpha 1$ -2 Fucosidase incubated for 20 hours at 37°C results in no detectable activity as determined by thin layer chromatography.

Glycosidase Activity (α 1-6 Galactosidase) - A 10 μ 1 reaction in Glyco Buffer 1 containing 1 nM of fluorescently-labeled α -Galactosidase substrate (Gal α 1-6Gl α 1-6Gl α 1-2Fru-AMC) and 100 units of α 1-2 Fucosidase incubated for 20 hours at 37°C results in no detectable activity as determined by thin layer chromatography.







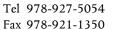

www.neb.com info@neb.com

New England Biolabs Product Specification

Assay Name/Specification (minimum release criteria)

Glycosidase Activity ($\alpha 1$ -6 Mannosidase) - A 10 μ l reaction in Glyco Buffer 1 containing 1 nM of fluorescently-labeled α -Mannosidase substrate (Man $\alpha 1$ -6Man $\alpha 1$ -6(Man $\alpha 1$ -3)Man-AMC) and 100 units of $\alpha 1$ -2 Fucosidase incubated for 20 hours at 37°C results in no detectable activity as determined by thin layer chromatography.

Protease Activity (SDS-PAGE) - A 20 μ l reaction in 1X Glyco Buffer 1 containing 24 μ g of a standard mixture of proteins and a minimum of 100 units of α 1-2 Fucosidase incubated for 20 hours at 37°C, results in no detectable degradation of the protein mixture as determined by SDS-PAGE with Coomassie Blue detection.

KruhKotum

Date 11 Mar 2016

Derek Robinson
Director of Quality Control





